


FOREHEAD THERMOMETER

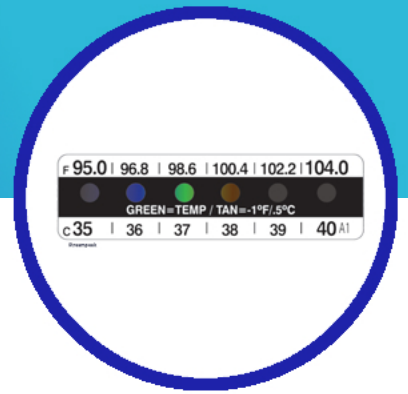

Product Description

During the outbreak of illness, early fever detection is one of the necessary precautions to take to prevent the further spread of illness and to ensure faster treatment.

Here in Stream Peak, we offer Forehead Thermometers to ensure the safety of you and your workers. Forehead Thermometers are technology-based liquid crystal temperature monitoring tool, which provides you with a simple yet cost-effective solution for your temperature monitoring needs.

This disposable thermometer allows you to monitor your temperature throughout the day by simply placing it on your forehead in indoors for approximately 15 seconds to obtain an accurate reading.


Industries Served

- Pharmaceutical

Product Usage

Forehead Thermometers should be used indoors and at room temperature. Direct sunlight and strong heat should be avoided when taking temperature. Temperature taking should not take place if the individual has been eating, drinking or exercising within the last 30 mins due to the possibility of temperature fluctuation.

1. Remove backing.
2. Hold thermometer on both ends and press it firmly on the dry forehead for about 15 seconds.
3. Read temperature while the thermometer is on the forehead
 - Green: Correct temperature
 - Blue: Add 1°F
 - Tan: Subtract 1°F

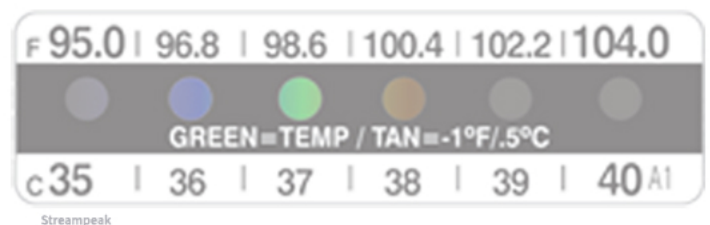


STREAM PEAK INTERNATIONAL PTE LTD
Tel: +65 6610 7921 | Fax: +65 6862 2520
Email: enquiry@streampeak.com.sg

Valuing Relationships, Protecting Goods
Address: 26 Tuas South Street 3, S(638026)
Website: www.streampeak.com.sg

Key Features

- Adhesive-backed paper
- Single or multi-use
- Fahrenheit & Celsius scale
- Made in the USA

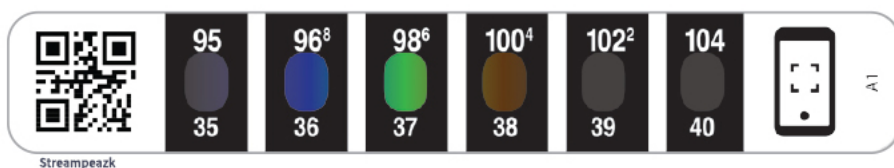

Dimensions

	P430 - Forehead Temperature Thermometer	C-8786 - Dual Scale Forehead Thermometer	C-8720 - Dual Scale Forehead Thermometer
Response Temperature	94°F, 96°F, 98°F, 100°F, 102°F, 104°F	35°C/95°F, 36°C/96°F, 37°C/98°F, 38°C/100°F, 39°C/102°F, 40°C/104°F	35°C/95°F, 36°C/96°F, 37°C/98°F, 38°C/100°F, 39°C/102°F, 40°C/104°F
Size	-	1.75" x .50"	2.825" x .44"
Application	Disposable; adhesive-backed	Disposable; adhesive-backed	Reusable; repositionable adhesive
Instruction for Use	-	-	Digitally available through QR code scan

Product Specifications

	P430 - Forehead Temperature Thermometer	C-8786 & C-8720 - Dual Scale Forehead Thermometer
Design Format	<ul style="list-style-type: none"> • Fahrenheit • Scale Hypoallergenic adhesive, pressure sensitive acrylate • Disposable 	<ul style="list-style-type: none"> • Fahrenheit and Celsius • Non-adhesive, medical grade adhesive, or hypoallergenic adhesive • Single use or repositionable
Range	94-104°F	Typically 95°F to 104°F and/or 35°C to 40°C
Accuracy	±1.0°F / 0.6°C	+/- 1.8°F / 1°C
Classification	FDA Class II Medical Device	FDA Class II Medical Device; Class I medical device in accordance with EU 2017/745. ISO 9001:2015 certified
Shelf Life	12 months	2 years
Configurations	-	Standard or Customizable
Packaging	Box of 100, case of 500	Packs of 100 or rolls of 5,000 with perforation

Related Products



DS System



Ice Gel Pack



FlashLink Mini
In-Transit Logger



Silica Gel



STREAM PEAK INTERNATIONAL PTE LTD
Tel: +65 6610 7921 | Fax: +65 6862 2520
Email: enquiry@streampeak.com.sg

Valuing Relationships, Protecting Goods
Address: 26 Tuas South Street 3, S(638026)
Website: www.streampeak.com.sg